


PROSTHETICS

STRAUMANN® METAL SCANBODIES

Straumann's new Metal Scanbodies feature durable titanium construction and careful design considerations to improve accuracy and support a precise digital workflow.

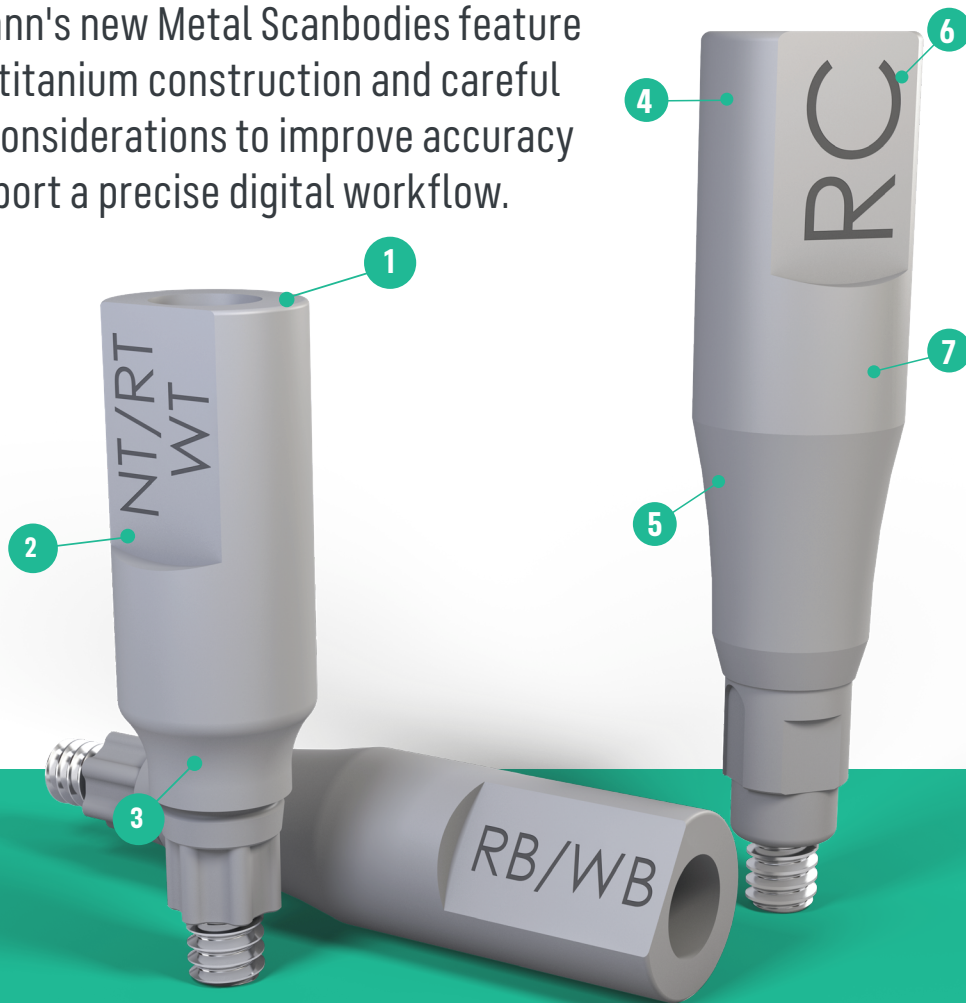

1. ROUNDED EDGES

No blackout in CAD software.

2. ONE FLAT SIDE

Helps define position in CAD software.

**3. HIGH AND CONTOURED
EMERGENCE PROFILE**

Reduces pressure on gingiva and optimizes visibility and scanability.

4. NARROW BODY

More space for adjacent teeth.

5. UNIQUE NON-REFLECTIVE COATING

Eliminates need for scan spray or powder.

6. LASER-MARKING

Easy to identify parts.

7. DURABLE TITANIUM


Radiopaque. Re-usable up to 100xs.

STRAUMANN® METAL SCANBODY

Quick reference


Implant Type	Platform	OLD PEEK Mono Scanbody	NEW Metal Scanbody
CrossFit® – Compatible to BLT Implant System	SC	025.0025	025.0077
	NC	025.2915	025.0078
	RC	025.4915	025.0079
CI RD – Compatible to Pure Ceramic Implant System	CI RD 4.1	032.041	032.130
SRA – Compatible to Screw-retained Abutments	SRA 3.5	025.0000	025.0080
	SRA 4.6	025.0001	025.0081
synOcta® – Compatible to Tissue Level Implant System	NNC	048.173	048.383
	RN	048.168	048.384
	WN	048.169	048.385
TorCFit™ BL – Compatible to BLC and BLX Implant Systems	RB/WB	065.0035	065.0103
TorCFit™ TL – Compatible to TLC and TLX Implant Systems	NT/RT/WT	036.3220	036.3230

CADCAM LIBRARY

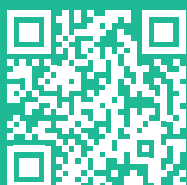


DOWNLOAD

CREATECH LIBRARY

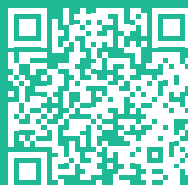

DOWNLOAD

LEARN MORE

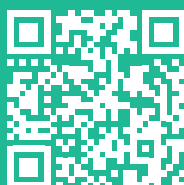

DOWNLOAD



Straumann Sweden
Straumann AB
Krokslätts Fabriker 45
SE-431 37 Mölndal
Sweden

+46 31 708 75 00
info.se@straumann.com
www.straumann.se

Straumann Denmark
Straumann Danmark ApS
Nygårds Plads 21, 1
DK-2605 Brøndby
Denmark

+45 46 16 06 66
info.dk@straumann.com
www.straumann.dk

Straumann Norway
Straumann AS
Nils Hansens vei 7
NO - 0667 Oslo
Norway

+47 23 35 44 88
info.no@straumann.com
www.straumann.no

Straumann Finland
Straumann Oy
Äyritie 12 A
01510 Vantaa
Finland

T +358 9 42 72 92 00
E info.fi@straumann.com
www.straumann.fi